Dental Solution

I. Introduction

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Dental Solution

I. Introduction

C2S thanks you for choosing the C2S pro. This machine and its patented system enables easier and risk-free **cleaning and upkeep of your removable dental appliances**.

The C2S pro has been designed and **developed by oral health specialists** to improve oral hygiene and comfort for denture wearers.

Developed for the **daily cleaning** of all types of removable dentures, its **innovative, fast and economic** cleaning system, which gently brushes the dentures with special cleaning needles, contributes to the **dentures' sustainability**. Use of individual, replaceable containers (every 6 months) ensures an **optimal hygiene** that improves comfort for the denture wearer and also prevents cross-contamination.

As with tooth brushing, denture maintenance with the C2S pro must be carried out daily for a good result. This machine is not intended to clean dentures with extensive tartar damage. If this is the case, then it is important to seek advice from a dental surgeon before using the C2S pro for the first time, to ensure that work begins on good denture bases.

IMPORTANT: THIS MACHINE MUST BE KEPT OUT OF THE REACH OF CHILDREN (UNDER 8 YEARS) AND MUST NOT BE USED BY THOSE WITH PACEMAKERS.

Contents of the C2S pro (CS-001) pack

1. A C2S pro machine
2. Two needle retaining magnets in the container (fig. 1)
3. An adapter connected to the machine
4. A user guide

Contents of the container pack: Whirlbowl (CS-002)

1. Two containers for the C2S pro (maximum use: 6 months)
2. Two pairs of tweezers
3. Two 40-g sachets of cleaning needles
4. A set of identification labels
5. A user guide

II. Product description


III. How to use

1/Turning on the machine

- Attach the power cable to the transformer.
- Connect the power cable to the machine and to a standard plug socket.
- Turn on the machine using the O/I switch behind the machine.
- The machine goes into sleep mode after 15 minutes of «non-use».
- To wake it up, press any of the buttons once.
- After each use, turn off the C2S using the O/I switch on the back of the machine.

2/Preparing the cleaning container

- Open the sachet of metal cleaning needles directly into the container ①. For better safety, position the needle retaining magnet under the container ②.
- Close the container ③. (Fig.1)

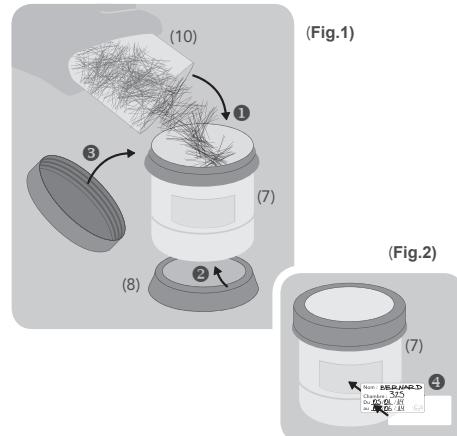4/ Programming and starting the cleaning cycle

- Press the «start/stop» button.
- You can stop the cleaning cycle at any time by pressing the «start/stop» button.
- During the cleaning, the motor will stop 10 seconds for every 3 minutes.
- Do not remove or open the container during the cleaning cycle.


Specific case: for the entire upper denture (upper jaw). We recommend turning the denture halfway through the cycle for a better clean. Example: for a 10-minute cleaning cycle, program the timer for 5 minutes for one side then 5 minutes for the other side.

- Use the «up/down» arrows to scroll to the correct time. The timer is pre-programmed for 10 minutes. The time reduces and increases in 5-minute increments up to 20 minutes.
- Timing selection depends on the amount of residue on the denture. Between 5 and 10 minutes for a standard upkeep clean (standard amount of residue). Between 10 and 15 minutes for quite a lot of residue. Between 15 and 20 minutes for a lot of residue.

Warning: if the tartar remains after a 20-minute clean, we recommend cleaning 3 times a day. If the tartar persists after a week of use, seek advice from a dental surgeon about removing the tartar and using the C2S pro on a clean base in order to prevent the formation of new residue.

5/Retrieving the denture at the end of a cleaning cycle

- Remove the container from the machine.
- Position the retaining magnet under the container ①.
- Open the container ②.
- Retrieve the denture using the tweezers supplied ③.
- Rinse the denture with copious amounts of water. (Fig.5)

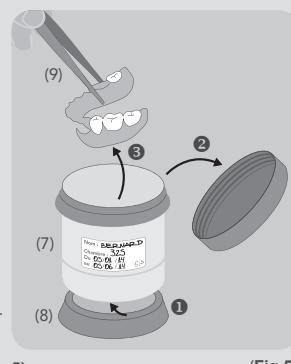
IMPORTANT: containers are strictly for individual use (this prevents cross-contamination) and must only be used on the C2S pro. Containers can be used for up to a maximum of 6 months. Should the needles fall out of their container (for example, if you forget to use the needle retaining magnet), it is prohibited to pick them up manually or to use the retaining magnet to gather them up. This could be dangerous and C2S is not liable for any problems that may be linked to such handling (bacteria on the needles, injuries etc.).

3/Positioning the container in the machine

- Position the container on the needle retaining magnet before opening it, in order to keep the needles at the bottom of the container ①.
- Open the container.
- Fill the container with tap water only, up to the level marking that is visible on the inside ②. For a better clean, we recommend using warm water.
- You can add a bicarbonate tablet. C2S recommends using this product to achieve a better clean. (Fig.3)

- Using the tweezers supplied, set the denture into the container ③.
- Close the container ④, remove the needle retaining magnet ⑤ and insert the container into the C2S pro. (Fig.4)

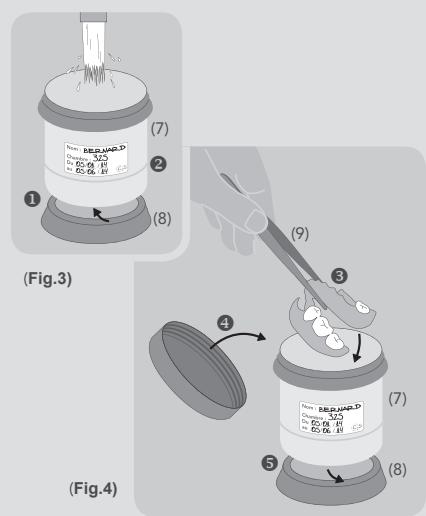
IMPORTANT: Thoroughly check that no metal needles are stuck in the denture!

- Dry the denture using a paper towel or place it in a glass of water to rinse it before putting it back in your mouth.
- Note:** after a long cleaning cycle, the water in the container may be quite hot. This is completely normal and is caused by the movement of the needles during the cleaning cycle.

6/Rinsing and maintaining the container after each use

- Keep the needle retaining magnet under the container.
- Empty the liquid in the container into a sink or waste water tank.
- Fill the container with water ① and use the retaining magnet to manipulate the needles in small circular motions ②. This will remove any residue remaining after the cleaning cycle. Repeat this step until the inside of the container is clean.
- Empty the remaining water from the container. (Fig.6)
- Close the container and remove the needle retaining magnet.
- Keep the container in a dry place that is protected from direct light.
- Replace the cleaning needles and container every 6 months.

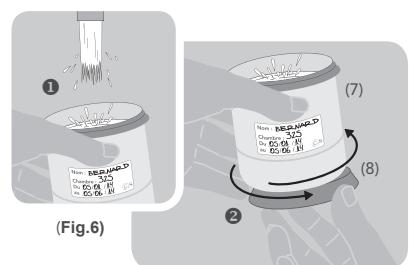Note: Use a damp cloth to clean the C2S pro. Do not use solvents or corrosive products.

- The used container and needles must be disposed of in specially-provided bins. The needles go in the needle bin and the container into the contaminated items bin.
- Do not use any needles that have been in another container or used for anything else.
- Do not re-use the same container without the needles.

For more information, go to www.praticima.fr.

WARNING

- The C2S must be kept out of the reach of children (under 8 years) and must not be used by those with pacemakers.
- Children must not play with this appliance.
- This appliance may be used by children who are 8 years or older and by those with reduced physical, sensory or mental abilities, or those lacking experience or knowledge, provided they are correctly supervised, or if the instructions related to the use of the appliance have been adequately given and if the risks involved have been understood.
- The C2S pro was designed for denture maintenance. As with tooth brushing, denture maintenance with the C2S pro must be carried out daily for a good result. This machine is not intended to clean dentures that are very damaged by tartar. If this is the case, then it is important to seek advice from a dental surgeon before using the C2S pro for the first time, to ensure that work begins on good tooth bases.
- Store the C2S pro in a dry place.
- Avoid prolonged exposure of the C2S pro to the sun or to any heat source.
- Use a damp cloth to clean the C2S pro. Do not use solvents or corrosive products.
- Children should not clean or maintain the appliance without supervision.
- After each use, turn off the C2S pro using the «O/I» switch behind the machine.
- The appliance must only be powered by the correct electrical supply.
- If the power cable is damaged, it must be replaced by the manufacturer, its After-Sales department or by similarly-qualified personnel, to prevent any risk.
- The transformer's flexible external cable cannot be replaced; should the cable become damaged, it must be disposed of and replaced with a transformer from the same model.
- The cleaning container must only be used with the C2S pro.
- Containers are strictly individual and must not be interchanged. If in doubt, changing the container is recommended in order to prevent any risk of cross-contamination.
- The metal needles and containers used for cleaning must never be separated for use.
- The container and needles must not be used (together or separately) for anything else other than as part of the C2S pro machine.
- The cleaning needles and container, as well as the C2S pro machine, should only be used for cleaning and maintaining dentures.
- The cleaning needles and container must be replaced every 6 months.
- Do not remove or open the container during the cleaning cycle.
- Should the needles fall out of their container (for example, if you forget to use the retaining magnet for the needles), it is prohibited to pick them up manually or to use the retaining magnet to gather them up. This could be dangerous and C2S is not liable for any problems that may be linked to such handling (bacteria on the needles, injuries etc.).
- Thoroughly check that none of the metal needles are stuck in the denture before putting it in your mouth.
- The used container and needles must be disposed of in specially-provided bins. The needles go in the needle bin and the container into the contaminated items bin.

Non-contractual photos of product improvements may appear.

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