

Neop | 600 x 400 single-column transfer cabinets

Supply of medicines or medical devices

- Robustness: excellent resistance, torsion-resistant steel frame, primer powder + polyester treated
- Possible binding for transport at all levels, on the external framework
- 100 % hygiene: designed with a minimum of retention zones, and resistant to disinfecting products of hospital environment
- Customization thanks to a wide range of drawers and accessories
- Ergonomic and handy: designed for medical staff comfort
- Unique design and harmonious colours



Robust



Hygienic



Ergonomic



Customizable

- Single-column transfer cabinet
- 2 heights: 160 or 180 cm
- Roller shutter door
- Code lock (with spare key) or badge lock
- 2 available inside panel models (slides) according to required internal fittings

Customizable structure

3 available colours

Magenta pink

Apple green

Sky blue



Standard configurations

Reference	Description
8YATR164H160_RM	NEOP "Basic" transfer cabinet 1 col. 600x400 H160 roller shutter door without lock - magenta pink
8YATR164H160_VP	NEOP "Basic" transfer cabinet 1 col. 600x400 H160 roller shutter door without lock - apple green
8YATR164H160_BC	NEOP "Basic" transfer cabinet 1 col. 600x400 H160 roller shutter door without lock - sky blue

Custom configurations

STEP 1 - Cabinet structure

Choice of cabinet structure according to inside panel models

	Cabinet height 160 cm*	Cabinet height 180 cm*
Standard inside panels 40 mm level	8Y164H160P4.TR	8Y164H180P4.TR
Optimized inside panels* 50 mm level	8Y164H160P5.TR	8Y164H180P5.TR

This reference includes:

- Structure
- Framework
- Wheels
- Inside panels


* for Modulo XL pill dispensers , 600 x 400 Modulo Classic and Eurojour pill dispensers

STEP 2 - Cabinet Customization

Selection of elements 1 to 3 is compulsory

		Cabinet height 160 cm*	Cabinet height 180 cm*	
1	Roller shutter door	8Y002_6016F3	8Y002_6018F3	
2	Code lock	8Y003CODE3	8Y003CODE3	
	Badge lock	8Y003BADGE3	8Y003BADGE3	
3	Colour	Magenta pink	8Y004_0016RM	8Y004_0018RM
		Apple green	8Y004_0016VP	8Y004_0018VP
		Sky blue	8Y004_0016BC	8Y004_0018BC
4	Wheels	8Y031	8Y031	

*for precise height, see next page

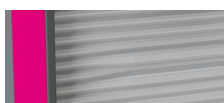

Roller shutter closure



Code lock



Badge lock

Magenta pink

Apple green

Sky blue

Accessories and internal fittings

Drawer, tow bar...

Please refer to product data sheets



Technical data



Materials

Structure	Complete welded one-piece structure, electro-galvanized steel, EPOXY anti-corrosion treated
Framework	Machine welded steel framework, powder primer coated + polyester treated
Slides	ABS (Acrylonitrile Butadiene Styrene)
Closing system	Roller shutter door: PP (Polypropylene), 2 mm thickness
Wheels	Semi-elastic rubber wheels - pressed steel bracket
Bumpers	4 stainless steel corners

Additional information

Wheels	160 mm diameter - 2 fixed wheels + 2 swivelling wheels with brake	
Empty weight (framework + structure)	H 160: 54 kg + 48 kg = 102 kg	H 180: 55 kg + 52 kg = 107 kg
Capacity	H 160 with standard inside panels (levels of 40 mm): 35 levels H 160 with optimized inside panels (levels of 50 mm): 27 levels	H 180 with standard inside panels (levels of 40 mm): 39 levels H 180 with optimized inside panels (levels of 50 mm): 30 levels
Locking system	Code lock / Badge lock	
Strip of colour	70 micron calendered PVC vinyl film, with acrylic adhesive	

Packaging

Packing	Individually on pallet strapped and film-wrapped, or on wheels, with proper protection	Dimensions L X D x H	H 160: 88 x 80 x 200 cm H 180: 88 x 80 x 215 cm
Gross weight (cabinet + packaging)	H 160: 114 kg H 180: 119 kg	Volume	H 160: 1,4 m ³ H 180: 1,5 m ³

Supply of medicines or medical devices

- Robustness: excellent resistance, torsion-resistant steel frame, primer powder + polyester treated
- Possible binding for transport at all levels, on the external framework
- 100 % hygiene: designed with a minimum of retention zones, and resistant to disinfecting products of hospital environment
- Customization thanks to a wide range of drawers and accessories
- Ergonomic and handy: designed for medical staff comfort
- Unique design and harmonious colours



Robust



Hygienic



Ergonomic



Customizable

- Double-column transfer cabinet
- Height: 160 cm
- Roller shutter door
- Code lock (with spare key) or badge lock
- 2 available inside panel models (slides) according to required internal fittings

Customizable structure

3 available colours

Magenta Pink

Apple green

Sky blue



Standard configurations - BASIC

Reference	Description
8YATR264H160_RM	NEOP "Basic" transfer cabinet 2 col. 600x400 H160 roller shutter door without lock - magenta pink
8YATR264H160_VP	NEOP "Basic" transfer cabinet 2 col. 600x400 H160 roller shutter door without lock - apple green
8YATR264H160_BC	NEOP "Basic" transfer cabinet 2 col. 600x400 H160 roller shutter door without lock - sky blue

Custom configurations

STEP 1 - Cabinet structure

Choice of cabinet structure according to inside panel models

	Cabinet height 160 cm*
Standard inside panels 40 mm level	8Y264H160P4.TR
Optimized inside panels* 50 mm level	8Y264H160P5.TR

* for Modulo XL pill dispensers , 600x400 Modulo Classic and Eurojour pill dispensers

This reference includes:

- **Structure**
- **Framework**
- **Wheels**
- **Inside panels**

STEP 2 - Cabinet Customization

Selection of elements 1 to 3 is compulsory

Closing system, locking system and colour strips must be double ordered for a double-column cabinet

		Cabinet height 160 cm*
1	Roller shutter door	8Y002_6016F3
2	Code lock	8Y003CODE3
	Badge lock	8Y003BADGE3
3	Colour	Magenta pink: 8Y004_0016RM Apple green: 8Y004_0016VP Sky blue: 8Y004_0016BC
4	Wheels	8Y031

*for precise height, see next page



Roller shutter closure



Code lock



Badge lock

Magenta pink

Apple green

Sky blue

Accessories and internal fittings

Drawer, tow bar...

Please refer to product data sheets



Technical data



Materials

Structure	Complete welded one-piece structure, electro-galvanized steel, EPOXY anti-corrosion treated
Framework	Machine welded steel framework, powder primer coated + polyester treated
Slides	ABS (Acrylonitrile Butadiene Styrene)
Closing system	Roller shutter door: PP (Polypropylene), 2 mm thickness
Wheels	Semi-elastic rubber wheels - pressed steel bracket
Bumpers	4 stainless steel corners

Additional Information

Wheels	160 mm diameter - 2 fixed wheels + 2 swivelling wheels with brake	
Empty weight (framework + structure)	64,4 kg + 94 kg = 158,4 kg	
Capacity	Cabinet with standard inside panels (40 mm levels) : 70 levels	Cabinet with optimized inside panels (50 mm levels) : 54 levels
Locking system	Code lock / Badge lock	
Strip of colour	70 micron calendered PVC vinyl film, with acrylic adhesive	

Packaging

Packing	Individually on pallet strapped and film-wrapped, or on wheels, with proper protection	Dimensions L X D x H	160 x 80 x 197 cm
Gross weight (cabinet + packaging)	170 kg	Volume	2,5 m ³