

# Neop | 400 x 400 emergency trolleys with telescopic drawers

## Material and products immediately available for first aid

- Robustness: excellent resistance , one piece structure made of aluminium
- 100% hygiene: designed with a minimum of retention zones, and resistant to disinfecting products of hospital environment
- Secure: locked with a numbered, single-use seal
- Ergonomic and handy: designed for medical staff comfort
- Modular: wide range of drawers and accessories
- Customisable: ward name identification on request



Robust



Hygienic



Ergonomic



Customizable

### Colour

Garance red

- Height 100 cm
- Ergonomic gripping area
- Door wing lock with 1 batch of 20 seals. Access to the content of the trolley breaking the seal by simple twist. After each use, the trolley is reequipped , then closed by a single use numbered seal .
- Drawers with full extension slides , with shock absorber and end stops, with perforated trays
- 5 drawer or 6 drawer version
- Frame for ISO trays adapted to most existing ISO trays
- Removable ISO tray for easy cleaning
- Cardiac massage board on support



## Standard configuration -PREMIUM


### EMERGENCY TROLLEY EQUIPPED WITH 5 telescopic drawers with compartmentalization

Treatment trolley 400 x 400 H100 equipped with:

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 1 telescopic drawer F400 H80 with full lip tray H50 equipped with a compartmentalization kit 16 boxes
- 3 telescopic drawers F400 H160 with openwork lip tray H100 equipped with a compartmentalization kit 9 boxes
- 1 telescopic drawer F400 H240 with openwork lip tray H100
- 1 support for cardiac massage board
- 1 cardiac massage board
- 1 monitoring support
- 1 IV pole with 2 hooks
- 1 lateral standard rail bottom right
- 1 oxygen tank holder

#### Reference

8Y144H100T.U2P

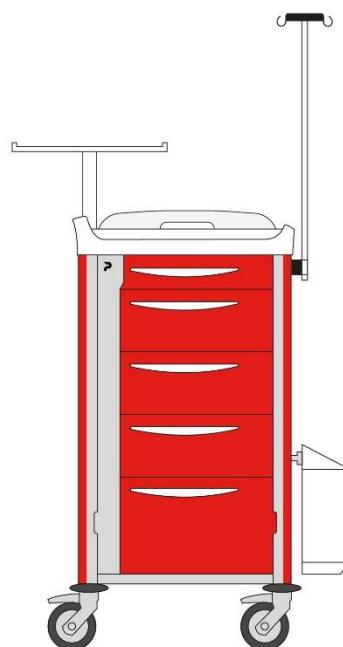


### EMERGENCY TROLLEY EQUIPPED WITH 5 telescopic drawers

Treatment trolley 400 x 400 H100 equipped with:

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 1 telescopic drawer F400 H80 with full lip tray H50
- 3 telescopic drawers F400 H160 with openwork lip tray H100
- 1 telescopic drawer F400 H240 with openwork lip tray H100
- 1 support for cardiac massage board
- 1 cardiac massage board
- 1 monitoring support
- 1 IV pole with 2 hooks
- 1 lateral standard rail bottom right
- 1 oxygen tank holder

#### Reference

8Y144H100T.U2

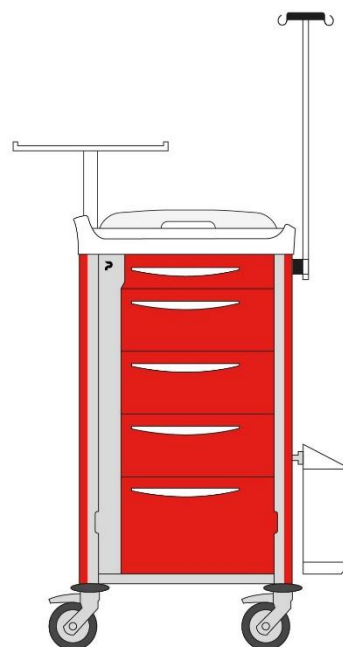

*Telescopic drawers fitted with lip trays can be divided with partitioning kits or combs, to be ordered separately.*


## EMERGENCY TROLLEY EQUIPPED WITH 6 telescopic drawers

Treatment trolley 400 x 400 H100 equipped with:

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 3 telescopic drawers F400 H80 with full lip tray H50
- 2 telescopic drawers F400 H160 with openwork lip tray H100
- 1 telescopic drawer F400 H240 with openwork lip tray H100
- 1 support for cardiac massage board
- 1 cardiac massage board
- 1 monitoring support
- 1 IV pole with 2 hooks
- 1 lateral standard rail bottom right
- 1 oxygen tank holder

### Reference

8Y144H100T.U4

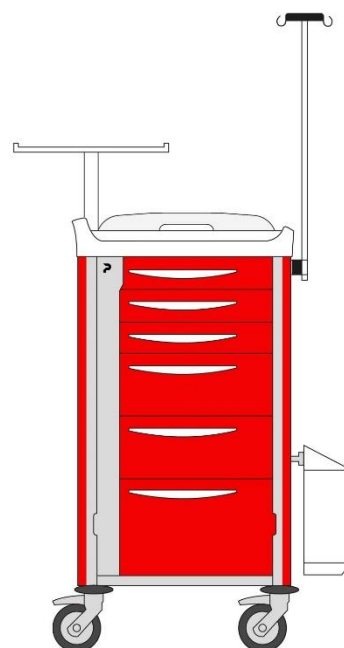

*Telescopic drawers fitted with lip trays can be divided with partitioning kits or combs, to be ordered separately.*


## BASIC EMERGENCY TROLLEY EQUIPPED WITH 5 telescopic drawers

Treatment trolley 400 x 400 H100 equipped with:

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 1 telescopic drawer F400 H80 with full lip tray H50
- 3 telescopic drawers F400 H160 with openwork lip tray H100
- 1 telescopic drawer F400 H240 with openwork lip tray H100
- 1 support for cardiac massage board
- 1 cardiac massage board

### Reference

8Y144H100T.U1

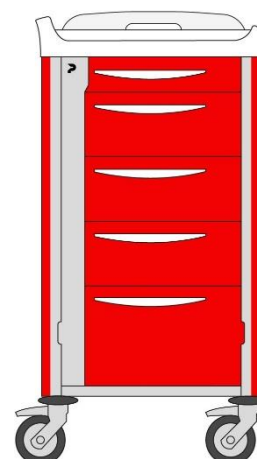

*Telescopic drawers fitted with lip trays can be divided with partitioning kits or combs, to be ordered separately.*


## BASIC EMERGENCY TROLLEY EQUIPPED WITH 6 telescopic drawers

Treatment trolley 400 x 400 H100 equipped with:

- 1 hinged closure with seal securing
- 1 batch of 20 seals
- 3 telescopic drawers F400 H80 with full lip tray H50
- 2 telescopic drawers F400 H160 with openwork lip tray H100
- 1 telescopic drawer F400 H240 with openwork lip tray H100
- 1 support for cardiac massage board
- 1 cardiac massage board

### Reference

8Y144H100T.U3

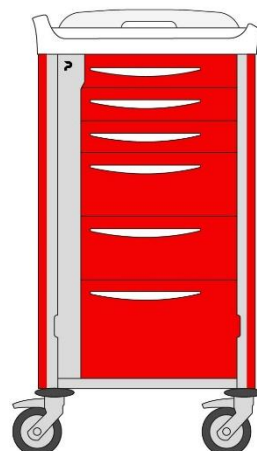

*Telescopic drawers fitted with lip trays can be divided with partitioning kits or combs, to be ordered separately.*

## Technical data

### Overall dimensions of the trolley

H. 1115 mm



The overall dimensions are measured with the wheels in the least bulky position.

### Useful dimensions of the worktop

L. 443 mm

D. 495 mm



Tablette supérieure en plastique thermoformé avec zone de préhension ergonomique

### Worktopheight

H. 1057 mm

## Materials

Structure	One-piece aluminium structure (20/10th) with continuous welding, treated with an epoxy anti-corrosion coating
Work surface	ABS / PMMA thermoplastic
Strip of colour	70 micron calendered PVC vinyl film with acrylic adhesive
Slides	Ballbearing slides made of galvanized steel with shock absorber
Closing system	Door wing: Aluminium with Epoxy painting
Wheels	Simple non marking rubber wheels - Polyamide cover (PA)
Bumpers	4 non marking round stop blocks made of EVA (Ethylene Vinyl Acetate)

## Additional information

Wheels	125 mm diameter – 2 breaking wheels + 2 swivelling wheels - wheel and wheel-pin double
Empty weight	35 kg
Locking system	Seals
Colors	70 micron calendered polymeric PVC film coated with acrylic adhesive

## Packaging

Packing	Individually on pallet strapped and film-wrapped, or on wheels, with proper protection
Gross Weight (trolley + packaging)	52 kg
Dimensions	L 0,80 x D 1,20 x H 1,45 m
Volume	1,15 m <sup>3</sup>