

Clean linen range

Nursing trolleys

Photos are not contractual

The range of nursing trolleys streamlines and facilitates the work of nursing assistants. The **ergonomic** design and **modular** concept of storage compartments make it easy to adapt the trolley to each ward's needs.

To limit the pushing effort required from users and prevent musculoskeletal disorders, our trolleys can be equipped with the **EasyMoove motorization system**. **Computerized management** is possible by adding a secured computer holder to the trolley.



Advantages

- ▶ Open or closed model with sliding doors
- ▶ 2 models: 15 to 20 patients and 20 to 30 patients
- ▶ Adjustable-height middle tray with removable dividers
- ▶ Several storage areas
- ▶ Materials resistant to corrosion and disinfecting products



Hygienic



Ergonomic

Clean linen range

Nursing trolleys

Mini nursing trolleys: 15 to 20 patients

Standard equipments

- ▶ 1 600 x 400 height 75 sliding drawer under upper shelf (only for open trolley)
- ▶ 1 adjustable and removable middle shelf
- ▶ 2 removable dividers for the middle shelf
- ▶ 4 standard rails
- ▶ 1 storage bin
- ▶ 2 large side baskets

| Description | Reference |
|---|-----------|
| Mini nursing trolley without doors | 8P500 |
| Mini nursing trolley without doors with upper divider | 8P500+ |
| Mini nursing trolley without doors with EasyMoove® motorization system | 8P500_EM |
| Mini nursing trolley with doors and code lock | 8P500_C |
| Mini nursing trolley with doors and key lock | 8P500_P |
| Mini nursing trolley with doors, key lock and upper divider | 8P500+_P |
| Mini nursing trolley with doors, key lock and EasyMoove® motorization system | 8P500_PEM |

Dimensions



Clean linen range

Nursing trolleys

Classic nursing trolleys: 20 to 30 patients

Standard equipments

- ▶ 1 600 x 400 height 75 sliding drawer under upper shelf (only for open trolley)
- ▶ 1 adjustable and removable middle shelf
- ▶ 3 removable dividers for the middle shelf
- ▶ 4 standard rails
- ▶ 1 storage bin
- ▶ 2 large side baskets

| Description | Reference |
|--|-----------|
| Classic nursing trolley without doors | 8P550 |
| Classic nursing trolley without doors with upper divider | 8P550+ |
| Classic nursing trolley without doors with EasyMoove® motorization system | 8P550_EM |
| Classic nursing trolley with doors and code lock | 8P550_C |
| Classic nursing trolley with doors and key lock | 8P550_P |
| Classic nursing trolley with doors, key lock and upper divider | 8P550+_P |
| Classic nursing trolley with doors, key lock and EasyMoove® motorization system | 8P550_PEM |

Dimensions



Clean linen range

Nursing trolleys

Technical data

Materials

| | |
|----------------------------|---|
| Structure | Structure: 40x20 steel tube - Shelves: 1.5mm thick welded steel sheet The entire assembly is treated with epoxy anti-corrosion coating |
| Wheels | Dual-band Polyamide (PA) wheels - Flexible Polyurethane (PUR) tread |
| Bumpers | 4 round bumpers made of non-marking EVA (Ethylene Vinyl Acetate) |
| Doors | Translucent PMMA, 6 mm thick |
| Removable shelf | Steel wire Ø 6 and Ø 4 treated with epoxy anti-corrosion coating |
| Adjustable dividers | 6 mm diameter steel wire with epoxy anti-corrosion treatment |
| Storage bin | ABS (Acrylonitrile Butadiene Styrene) |
| Drawer | ABS injection-molded drawer 600 x 400 |
| Basket | Steel coated with epoxy polyester RAL7024 |

Additional information

| | | | | |
|-----------------------|---|--------------------------------|----------------------------------|------------------------------------|
| Wheels | 4 swivelling wheels Ø 125 mm (1 with brake and 1 with directional lock) | | | |
| Pre-equipment | To be struck on the plateau | | | |
| Empty weight | 8P500: 38 kg 8P550: 51 kg | 8P500+: 41 kg 8P550+: 57 kg | 8P500_P: 51 kg 8P550_P: 65 kg | 8P500+_P: 54 kg 8P550+_P: 68 kg |
| Load capacity | Top shelf: 60 kg - Removable shelf: 30 kg - Lower shelf: 60 kg | | | |
| Locking system | For trolleys with doors: key lock or code lock (with spare key) | | | |

Packaging

| | | |
|---------------------|--|--------------|
| Packing | Individually on pallet + cardboard packaging | |
| Gross weight | 8P500: 50 kg | 8P550: 66 kg |
| Dimensions | L. 0,80 x D. 1,20 x H. 1,30 m | |
| Volume | 1,3 m ³ | |

Colors

| <i>Description</i> | | <i>Reference</i> |
|---------------------------|---|------------------|
| Pop blue color |  | 8P504_BP |
| Magenta pink color |  | 8P504_RM |
| Blueberry color |  | 8P504_MY |
| Orange color |  | 8P504_OR |
| Apple green color |  | 8P504_VP |

Clean linen range

Nursing trolleys

Motorization



All our nursing trolleys can be equipped with the motorized wheel EasyMoove®. This innovation enables to **reduce physical efforts of the users**. A simple act with the finger on the handle allows the system to operate. Staff's work is therefore with utmost convenience.

| Description | Reference |
|--|-------------|
| EasyMoove® motorization system for Mini nursing trolley (to be ordered at the same time as the trolley) | 9R002_L02V2 |
| EasyMoove® motorization system for Classic nursing trolley (to be ordered at the same time as the trolley) | 9R002_L03V2 |



IT supports

| Description | Reference |
|---|-----------|
| 1 Sheet metal laptop holder | 8Y009CP4 |
| 2 PMMA laptop holder with cover | 8Y009CP2 |
| 3 PMMA laptop holder without cover | 8Y009CP1 |
| 4 Touch tablet holder | 8Y009CP3 |
| 5 Vesa screen holder | 8Y009CP5 |
| 6 Computer arm extension | QUI2F110 |
| Compact mouse holder for compact laptop holder | 28502 |
| Cable guide | 8E105 |
| 10 cable ties | 8E109 |

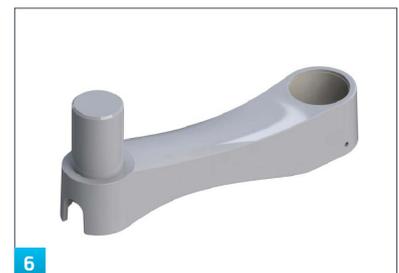
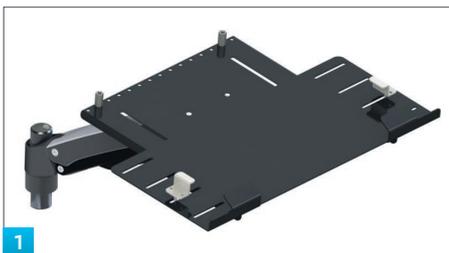


Clean linen range

Nursing trolleys

Accessories for standard rails

| | Description | Overall dimensions (in mm) | Reference |
|----|--|----------------------------|-----------|
| | Side basket small model | L. 355 x P. 218 x H. 238 | 28062N.GR |
| 7 | Side basket big model | L. 405 x P. 218 x H. 238 | 28061N.GR |
| 8 | Storage bin | L. 455 x P. 135 x H. 110 | 28063N.GR |
| 9 | Chrome wire briefcase | | BF20250 |
| 10 | Document sorter - 3 compartments | L. 257 x P. 130 x H. 445 | 28011N |
| 11 | Vertical support for allocation sheets | | AN09403 |
| 12 | Waste bin holder | L. 173 x P. 265 x H. 302 | 28080.GR |
| 13 | Pedal waste bin 25 L | L. 235 x P. 190 x H. 400 | 28074N |
| 14 | Standard bracket for gloves | L. 135 x P. 80 x H. 300 | 28001N.GR |
| 15 | Removable roll support | | 28059N.GR |

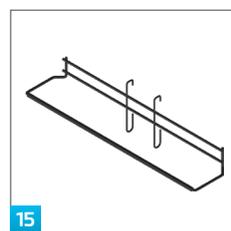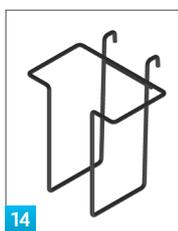

Other equipment

| | Description | Reference |
|----|---|-----------|
| 16 | Handle | 8P502 |
| 17 | Removable 3-compartment upper divider | 8P515 |
| 18 | Bottle holder - small model | 28038 |
| | Removable shelf for Mini nursing trolley | ACC8P505 |
| | Removable shelf for Classic nursing trolley | ACC8P554 |
| 19 | Fixed bag holder with manually opening lid | 8P525 |
| 20 | Removable bag holder with lid and opening pedal | 22001_S |

